
For a full explanation of all features and instructions, please refer to the User's Guide (available for download from support.hubbleconnected.com).

1. Setting up your Smart Bluetooth Baby Scale

A. Install battery

The battery operates with 4 AAA alkaline batteries type LR03 (1.5V).


- Turn the weighing scale bottom up to locate the battery compartment, open the battery door to install the batteries provided as shown in figure.
- Install four (4) new "AAA" batteries (LR03) into the battery compartment according to the polarity markings inside.
- Replace the battery cover.
- Press the Power ON/OFF button to turn on the weighing scale.

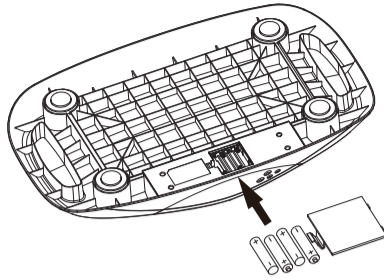Note:

- Batteries are to be installed with the correct polarity.
- Do not mix batteries of different types or mix old and new batteries, replace all 4 batteries when empty with new batteries.
- If the scale is to be stored unused for a long period of time, the batteries should be removed.
- The battery contacts are not to be short-circuited.

B. Set up HubbleClub by Hubble Connect App to track your baby's growth

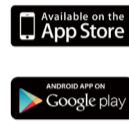
- Scan the QR code with your smart device and download HubbleClub by Hubble Connected App from the App Store for iOS devices or from the Google Play™ Store for Android™ devices.
- Install HubbleClub by Hubble Connected App on your device.

- Open the HubbleClub by Hubble Connected App on your compatible smartphone or tablet.
- Turn on Bluetooth on your compatible smartphone or tablet.
- Follow the in-app instructions to create your Hubble account and connect to your device.

Note:

Please take note of the following minimum system requirements: Smartphones/Tablets: iOS 10.0, Android™ 7.0, Bluetooth 4.2.

2. Overview of the control panel

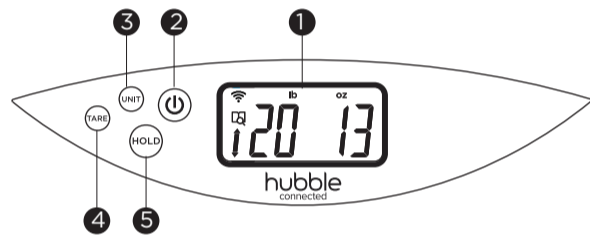1. LCD Display

- Bluetooth Signal
- Memory symbol
- Flashes when the memory is full.
- Weight Change indicator
- The symbol appears together with the Memory symbol to show weight value change with Up arrow (weight gain) or Down arrow (weight loss) indications.

2. Power ON/OFF button

- Select the unit of measurement: kilograms (kg) or pounds (lb).

3. Using the baby scale

3.1 Turn On/Off the baby scale

Press the Power ON/OFF button to turn On/Off the weighing scale.

3.2 Select Weighing Unit

Press the UNIT button on the control panel to select unit of measurement, kilograms (kg) or pounds (lb).

3.3 Tare Weight

The function is to reset the scale's display to zero. When you measure your baby together with the soft pad, you don't want to include the soft pad weight in the reading.

- To tare the soft pad, place it on scale, wait for a stable reading.
- Press the TARE button to reset the display to zero, the digit "0" is displayed.


3.4 Baby Weighing

- Press Power ON/OFF button to turn on the baby scale, when the digits "0" are displayed, the scale is ready for weighing.
- Place your baby on the weighing scale until the weight reading steadily appears on the display.
- Press the HOLD button to lock and store the data. The data is transmitted to HubbleClub app when smart device and scale are in sync.


3.5 Weight Tracker

The scale tracks the weight change between current and last measurements.


The display shows the current weight and the weight difference (between the current weight and the last weight) alternately. Each weight displays for 2s. The scale will turn off automatically after 8s. (Note: In the first use after installing the batteries, the display will only show the current weight, then turn off automatically after 8s.)

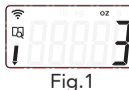- The symbol ↓ indicates it is weight lose against last reading.

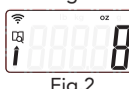- The symbol ↑ indicates it is weight gain against last reading.

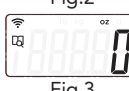- No weight change was detected (Fig.3).

Note:

In the first use after installing the batteries, the display will only show the current weight. The weight data will not be kept whenever the batteries are replaced.

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